

GIO perfect Ultrasonic Skin Scrubber with EMS Therapy



USER MANUAL

To avoid harm caused by misuse, please read the instructions carefully and use this product strictly in accordance with the manual.

PRODUCT SPECIFICATION

Product name: GIOperfect Ultrasonic Skin Scrubber with EMS Therapy
Product Model: GIO-730
Dimensions: 130 × 70 × 17 mm
Weight: 75 g
Battery Voltage: 3.7 V
Rated Power: 1.5 W
Battery Capacity: 250 mAh
Blade Frequency: 26,000 vibrations/second
Charging Time: Approx. 1–1.5 hours
Working Time: Approx. 0.8 hour

PACKAGE CONTENTS

- 1 × Device
- 1 × Charging Cable
- 1 × User Manual

SAFETY INFORMATION

- The warning labels indicate potential danger. Failure to follow instructions may result in injury or property damage.
- Stop using immediately if you experience pain, discomfort, or abnormal reactions.
- Always turn off and unplug the device before cleaning or maintenance.
- Keep out of reach of children.
- Do not use the device near flammable substances (e.g., alcohol, gasoline, acetone, sprays).
- Do not immerse the product in water or rinse it directly under running water.
- If the device is dropped or damaged, do not continue to use it.
- Do not disassemble or attempt to repair the device yourself. Contact the original point of purchase for assistance.
- Only use the charging cable provided and ensure the power supply meets the device's specifications.
- Clean the device after each use to prevent residue build-up.
- Do not use the device during a bath or shower, or near bathtubs, basins, or other water containers.

Not suitable for use by:


- Pregnant or menstruating individuals, or those with hormone-related skin issues.
- Individuals with allergies, skin diseases, dermatitis, or severe sunburn.
- Those with open wounds, eczema, or swelling.
- Individuals recovering from surgery.
- People with physical, sensory, or mental impairments, unless supervised by a guardian.

PRODUCT OVERVIEW

Proper skin cleansing is essential for product absorption. Without thorough cleaning, skincare products may not penetrate the skin effectively, leading to dullness, acne, sagging, dark spots, and other skin issues.

This device utilizes high-frequency ultrasonic vibration and ion technology to deeply cleanse the skin and enhance absorption of skincare products.

USB-C charging cable

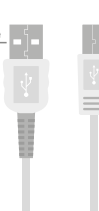

Ion Contact Bar



Dust cover

Power / Mode Button




HOW TO USE

PREPARATION BEFORE USE

1. Charge the device for 2 hours before first use.
2. Press and hold the power button for 2 seconds to turn on/off. Short press to switch modes.
3. Clean the blade with a soft cloth or cotton pad.
4. Apply a cleanser or clean water to the skin.
5. When the battery is low, the device will emit two beeps and shut down.
6. During charging, indicator lights will flash one by one. When fully charged, all lights will remain on.


1. Ultrasonic Cleansing

- Wash your face and apply facial cleanser.
- Press the first mode to activate vibration cleansing.
- Hold the device and move the blade slowly at a 45° angle against the direction of pores.

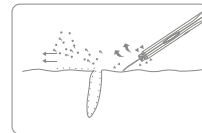

2. EMS Detox

- Cleanse and steam your face for 5 minutes (or apply a warm towel).
- Apply blackhead removal solution for 10 minutes.
- Press the fourth mode to activate positive ion export.
- Move the blade slowly downward, against the direction of pores.
- Use for 10 minutes, then rinse face and clean the device.

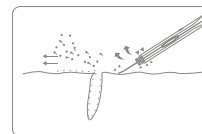

3. Nutrient Infusion

- Cleanse and steam the face for 5 minutes.
- Apply a nutrient-rich essence.


Mode 3 (Fast infusion):

- Press the third mode (negative ion import).
- Move the blade slowly from bottom to top, along the direction of pores. Use for 5 minutes.



Mode 2 (Slow patting infusion):


- Press the second mode to activate patting massage.
- Again move slowly from bottom to top. Use for 5 minutes.
- After treatment, rinse the face and steam with a cold spray or apply a cold towel for 10 minutes to tighten pores.

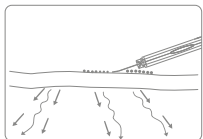

CLEANING PROCEDURES



Face cleaning and care

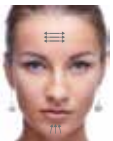

Forehead and jawline care



Nose and T-zone care

1. Focus on key areas: forehead, jawline, and T-zone (including the nose bridge).
2. Always hold the ion contact bar on the side of the device to activate ion functions.
3. Keep the blade at a 45° angle and move it slowly and gently across the skin — do not leave it stationary.
4. Do not use on inflamed or irritated skin.
5. Recommended usage: 2–3 times per week.
6. Always clean the device and blade before and after use.

CLEANING & MAINTENANCE

- Store the device in a dry place. Do not place it in dishwashers or microwaves.
- Turn off and unplug before cleaning.
- Use only soft cloths dampened with cold or warm water — avoid sharp tools, corrosive agents, or volatile liquids.
- Do not rinse the device directly under water.
- If unused for more than 3 months, fully charge the device before storage to prevent battery damage.

TROUBLESHOOTING

Issue	Possible cause	Solution
Shovel doesn't vibrate	Shovel is deformed	Gently hold the shovel with both hands and adjust it back into shape
	Internal wiring is disconnected	Contact after-sales customer service
	Excessive pressure applied during use	Use light pressure — gently touch the skin with the shovel
Device won't power on	Low battery	Charge the device for at least 2 hours before use
	Battery is damaged due to long period of disuse	Contact after-sales customer service
	Water has entered the device	Contact after-sales customer service
Ion function failure	Moisture has entered the power button	Contact after-sales customer service
	User's hands are not in contact with the side electrode ("contact bar")	Ensure your hand is touching the lateral electrode during use
	Internal circuit damaged by water ingress	Contact after-sales customer service

- 🇬🇧 Scan the QR code to download the manual
- 🇨🇵 Zeskanuj kod QR, aby pobrać instrukcję
- 🇩🇪 Scannen Sie den QR-Code, um das Handbuch herunterzuladen
- 🇫🇷 Scannez le code QR pour télécharger le manuel
- 🇪🇸 Escanee el código QR para descargar el manual
- 🇵🇹 Digitalize o código QR para descarregar o manual
- 🇮🇹 Scansiona il codice QR per scaricare il manuale
- 🇷🇺 Scan de QR-code om de handleiding te downloaden
- 🇺🇾 Скануйте QR-код, щоб завантажити інструкцію
- 🇨🇿 Naskenujte QR kód a stáhněte si návod k použití
- 🇸🇰 Naskenujte QR kód a stiahnite si príručku
- 🇷🇴 Scanează codul QR pentru a descărca manualul
- 🇵🇱 Nuskaitykite QR kodą, kad atsisiųstumėte instrukciją
- 🇱🇻 Nolasiet QR kodu, lai lejupielādētu rokasgrāmatu
- 🇩🇪 Lae juhend alla, skannides QR-koodi
- 🇭🇺 A kézikönyv letöltéséhez olvassa be a QR-kódot.
- 🇵🇪 لوليدال ليزناتل قعي رسال قباچت سالا اذر حسب



innoGLO

Manufactured in China for
INNOGIO GLOBAL sp. z o.o.
Ostródzka street 74H, 03-289 Warsaw, Poland
www.innogio.com

@innogio.official @innogio_com

