

EN

GIO-731

innoGIO

GIO perfect Skin Scrubber


Ultrasonic face scrubber



USER MANUAL

To avoid harm caused by misuse, please read the instructions carefully and use this product strictly in accordance with the manual.

PRODUCT OVERVIEW

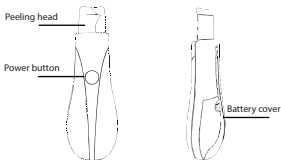

SAFETY & CARE INSTRUCTIONS

Before Use

- Insert the battery (2xAAA, not included) and close the battery cover securely.
- Make sure the device head is clean and dry before using.

During Use

- Use only on moist skin — never on dry skin.
- Recommended usage: 1–2 times per week, for up to 10 minutes per session.
- Do not use on broken, inflamed, irritated, or sunburned skin.
- Discontinue use immediately if you experience pain or discomfort.
- Avoid using metal objects (e.g. jewelry, keys) near the device.
- This device is not waterproof — do not use in the shower or bath.

After Use

- Clean the device head with a tissue or soft cloth immediately after each use.
- Ensure no lotion, gel, or moisture remains on the device before storing.
- Store in a cool, dry place away from direct sunlight, heat sources, and humidity.
- Keep the device safely stored and out of reach of children.

HOW TO USE

Always moisturize your skin before starting!

Step-by-Step:

- Remove makeup and fully hydrate your face with a spray or mist.
- Warm the skin using a hot towel for improved effectiveness.
- Hold the device ensuring your hand does not contact the metal plate.
- Tilt the device to a 30°–45° angle.
- Move gently in the direction shown (see below) — always glide, do not press.
- After use, power off and clean the device.
- Apply toner or serum to finish your skincare.

How to move the peeling head:

<u>Facial Area</u>	<u>Direction of Use</u>
Chin	Move upward
Nose	Glide bottom to top
Cheeks	Move outward and upward from center
Mouth area	Glide towards center of face

Do Not Use On:

- Around or on the eyes
- Eyelids
- Infected, injured, or ultrasensitive skin



-  Scan the QR code to download the manual
-  Zeskanuj kod QR, aby pobrać instrukcję
-  Scannen Sie den QR-Code, um das Handbuch herunterzuladen
-  Scannez le code QR pour télécharger le manuel
-  Escanee el código QR para descargar el manual
-  Digitalize o código QR para descarregar o manual
-  Scansiona il codice QR per scaricare il manuale
-  Scan de QR-code om de handleiding te downloaden
-  Скануйте QR-код, щоб завантажити інструкцію
-  Naskenujte QR kód a stáhněte si návod k použití
-  Naskenujte QR kód a stiahnite si príručku
-  Scanează codul QR pentru a descărca manualul
-  Nuskaitykite QR kodą, kad atsisiųstumėte instrukciją
-  Nolasiet QR kodu, lai lejupielādētu rokasgrāmatu
-  Lae juhend alla, skannides QR-koodi
-  A kézikönyv letöltéséhez olvassa be a QR-kódot.
-  ليلا هذا ليؤنزل فيجيسر الكايفاتسالة زمر جسرما



innoGIO

Manufactured in China for **INNIOGIO GLOBAL** sp. z o.o.
 Ostródzka street 74H, 03-289 Warsaw, Poland
www.innogio.com



@innoGIO.official



@innogio_com