


EN

GIO-750

innoGIO



GIO perfect Neck Brace

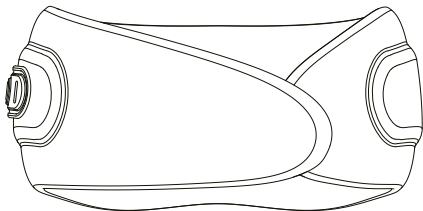

USER MANUAL

To avoid harm caused by misuse, please read the instructions carefully and use this product strictly in accordance with the manual.

PRODUCT OVERVIEW

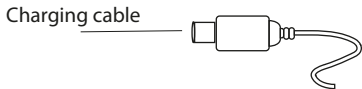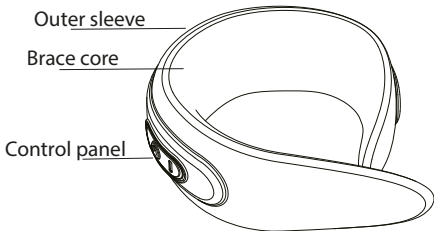

PACKAGE CONTENTS

- Main unit
- USB-C charging cable
- User manual

OPERATING INSTRUCTIONS

1. Press and hold the power button for 1.5 seconds to turn on (green light indicates power).
2. Short press to cycle through 3 modes (green, blue, red).
3. Double-click to activate/deactivate the heating function (55°C; off by default).
4. Heating and vibration functions can operate simultaneously.
5. Device auto-shuts off after 15 minutes. Restart manually if needed.
6. Press and hold for 1.5 seconds to turn off during use.
7. Recharge when the battery is low.
8. Use only the original charging cable.
9. Do not use the device while charging.
10. A full charge takes 3–4 hours (green light indicates completion). Unplug before use.
11. Accessories are designed exclusively for this device.

TROUBLESHOOTING

Issue	Possible Cause	Solution
Device won't start	Low battery	Charge and restart
Stops during use	Timer ended or battery low	Recharge and restart

IMPORTANT SAFETY INFORMATION

- This product includes a heating function — use with caution, especially if you are heat-sensitive.
- Charge the battery at least once every 1–2 months to maintain optimal performance.
- Use only a certified 5V 1A charger. Higher voltages may damage the device.
- Ensure the input voltage matches the product specifications.
- Do not use non-original accessories or substitutes.
- Do not charge the device during voltage fluctuations.
- Do not charge or use the product in bathrooms or humid environments.
- Avoid contact with the power plug using wet hands.
- Avoid prolonged exposure of the device to direct sunlight or high temperatures.
- Children or persons with reduced physical, sensory, or mental capabilities must use the device under supervision.
- Do not dispose of the built-in lithium battery with household waste. Follow local environmental regulations for proper disposal.



PRODUCT PRECAUTIONS

Do not use if:

- It fails to start.
- The battery depletes during standby.
- The session ends due to the timer or battery drain.

Consult a doctor before use if you:

- Are receiving medical treatment or feel unwell.
- Have an acute or chronic illness.
- Are under the influence of alcohol.
- Are a minor.

Battery Charging and Disposal

- Charge the battery at least once every 1–2 months to maintain performance.
- Use a certified 5V 1A charger. Higher voltages may damage the device.
- Do not dispose of the battery with household waste.

SPECIFICATION

Product Name: GIOperfect Neck Brace

Model: GIO-750

Dimensions: 180 × 210 × 90 mm

Power Consumption: 8W

Charging Voltage: 5V

Charging Current: 0.5A

Battery Capacity: 1200mAh

Charging Time: 3–4 hours

-  Scan the QR code to download the manual
-  Zeskanuj kod QR, aby pobrać instrukcję
-  Scannen Sie den QR-Code, um das Handbuch herunterzuladen
-  Scannez le code QR pour télécharger le manuel
-  Escanee el código QR para descargar el manual
-  Digitalize o código QR para descarregar o manual
-  Scansiona il codice QR per scaricare il manuale
-  Scan de QR-code om de handleiding te downloaden
-  Скануйте QR-код, щоб завантажити інструкцію
-  Naskenujte QR kód a stáhněte si návod k použití
-  Naskenujte QR kód a stiahnite si príručku
-  Scanează codul QR pentru a descărca manualul
-  Nuskaitykite QR kodą, kad atsisiųstumėte instrukciją
-  Nolasiet QR kodu, lai lejupielādētu rokasgrāmatu
-  Lae juhend alla, skannides QR-koodi
-  A kézikönyv letöltéséhez olvassa be a QR-kódot.
-  لي ليدل لي زنتل لة يرسلا لة باجات سالا زمر ح س ا



innogIO

Manufactured in China for **INNOGIO GLOBAL** sp. z o.o.
Ostródzka street 74H, 03-289 Warsaw, Poland
www.innogio.com



@innogIO.official



@innogio_com