

EN

GIO-780

innoGIO

# GIO perfect Mask


LED Facial Mask



**USER MANUAL**

To avoid harm caused by misuse, please read the instructions carefully and use this product strictly in accordance with the manual.


## GETTING STARTED

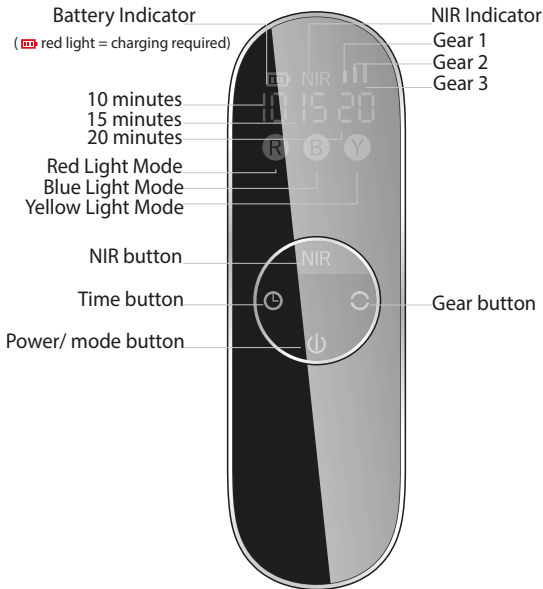

**Note:** When a controller symbol is displayed, the function is active.

### BUTTON FUNCTIONS

- Power/Mode button: Power ON/OFF and switch between light modes
- Time button: Select treatment duration (10, 15, or 20 minutes)
- Gear button: Adjust intensity (1–3 levels)
- NIR button: Activate/deactivate near-infrared lights

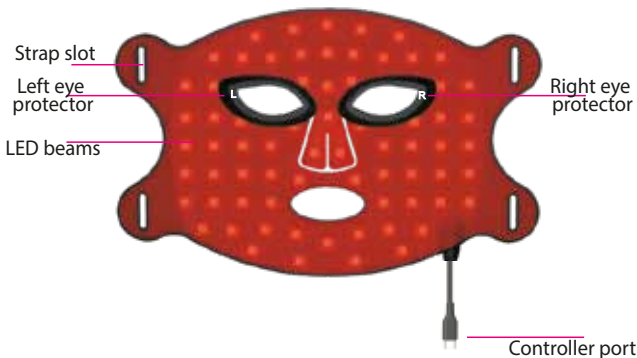

## ATTENTION

1. The device cannot be operated while charging.
2. Always clean the device before use with the power switched off. Use only recommended cleaning products (clean with a damp cloth or moist tissue, or use a saline solution (0.9% NaCl) to wipe parts that come into contact with the skin).
3. Store the LED Therapy Face Mask flat, in a dry place, at room temperature, away from direct sunlight.
4. Temporary redness may occur after treatment; it usually subsides within 1–2 hours.
5. After 12 treatments, users typically report visible improvement in skin tone and texture.
6. Fully charge the device with the supplied USB cable before first use. The controller light flashes during charging.

## INTENDED USE

- **Red light** (660 nm): Anti-aging, stimulates collagen, promotes relaxation, improves elasticity.
- **Blue light** (410 nm): Eliminates acne-causing bacteria, reduces inflammation and scars, increases oxygenation, brightens skin tone.
- **Yellow light** (590 nm): Balances skin texture, stimulates red blood cells, reduces redness/swelling/fine lines, fades acne marks, stabilizes skin condition.
- **NIR light** (850 nm): Repairs skin, accelerates metabolism, improves blood circulation, increases hydration.

## TREATMENT RECOMMENDATIONS

1. For best results, use 3–5 times per week for at least 4 weeks.
2. Maintenance treatments may be continued as needed, depending on skin condition and lifestyle.


## SENSITIZATION TEST

Before your first treatment, conduct a skin sensitivity test:

1. Place the LED Therapy Face Mask on your inner arm.
2. Press the Power Key to start a full 10-minute treatment.
3. If your skin feels hot or uncomfortable, stop immediately and remove the device.
4. Wait at least 6 hours. If redness or itching occurs, the device may not be suitable for you.


## HOW TO USE

### STEP 1

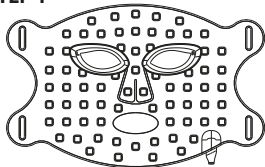

Remove the protective film

### STEP 2

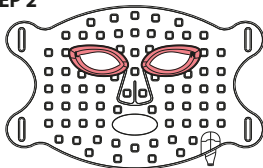


Install left and right eye protectors

### STEP 3



Install left and right straps

### STEP 4

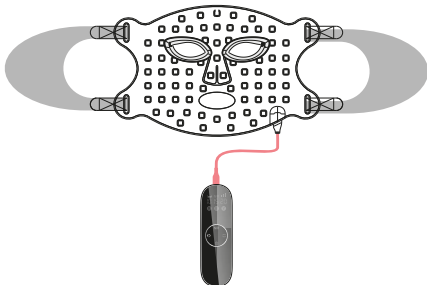

Type C plug-in controller

## OPERATING INSTRUCTIONS

### **Establish Skin Tolerance**

- Start with Gear 1–2 intensity.
- Gradually increase as your skin adapts.

### **Adjust Duration**

- Choose 10, 15, or 20 minutes with the Time button.
- Do not increase time or frequency excessively.

### **Recommended Usage**

- 3–5 times per week, 10–20 minutes per session.
- For Blue Light Therapy, do not exceed 10 minutes.

### **Carrying out a Treatment**

- Cleanse and exfoliate your skin.
- Place the mask, secure straps, and press the Power button.
- It starts in Red Light mode. Switch modes with the Power button, adjust intensity with the Gear button, and set time with the Time button.
- The mask turns off automatically.
- After use, remove the mask, apply skincare, and use sunscreen.
- Wipe clean with a damp cloth or water-based wipe. Do not clean while charging.

## PACKAGE CONTENTS

LED Mask ×1  
Controller ×1  
Elastic Straps ×2  
User Manual ×1  
Charging Cable ×1  
Storage Bag ×1  
Eye Protectors ×2

## PRECAUTIONS

- Use only for the conditions specified in this manual.
- Do not use after recent injury, surgery, or cosmetic procedures (e.g., fillers, Botox, microneedling, laser, chemical peel) until fully healed.
- Consult a doctor before use if pregnant, or if you have cancer or suspect you may have it.
- Do not treat red or swollen skin caused by previous treatment.
- Discontinue use if redness or swelling persists for more than 2 hours.
- Do not use if you have extremely sensitive skin, open wounds, infections, or herpes outbreaks.
- Do not use on painful, raised, crusting, or bleeding areas.
- Do not apply retinol before Red Light therapy (use only after).
- Some users may experience temporary redness, dryness, or irritation, especially with Blue Light.
- Protect skin with sunscreen after treatment, as sensitivity to sunlight may increase.
- Do not overuse Blue Light therapy; excessive use may cause pigmentation.
- Do not look directly at the LEDs. Use the supplied eye protectors.
- For adults over 18 years old; not recommended for children under 12.
- Do not modify the device.

## WARNINGS

### **Before Use:**

- Always conduct a sensitivity test before first treatment.
- Use only as specified; untested uses may pose unknown risks.

### **Device Safety:**

- Stop immediately if skin feels uncomfortably hot.
- Remove creams, serums, lotions, and makeup before treatment.
- Keep away from water, heat, and flammable areas.
- Do not modify or disassemble the device.
- Do not use if the device, battery, or charger is damaged.

## Charging Safety:

- Charge indoors, in a dry, ventilated area, with a USB-C charger.
- Do not charge near water sources or in explosive environments.
- Do not overcharge; disconnect after full charge (approx. 1 hour).
- Only use supplied accessories.
- Store in a cool, dry place.
- Charge between 0°C–40°C (32°F–104°F).

## Battery Safety:

- Do not crush, drop, or expose the battery to damage.
- Avoid short-circuiting the battery.
- If battery fluid contacts skin or eyes, rinse immediately and seek medical attention.

## Disposal

- Do not dispose of the device in household waste.
- The controller contains a lithium-ion battery. At end of life, return the unit to the manufacturer or follow local recycling and disposal guidelines for electronic devices.

## SPECIFICATION

Product name: GIOperfect Mask

Product ID: GIO-780

LED Colors: Red, Blue, Yellow, Infrared

Number of LEDs: 70 pcs (280 beams of light)

Battery: 2000 mAh

Voltage: 3.7V

Power: Max. 6.5W

Controller Size: 130 × 40 × 26.5 mm

Mask Size: 299 × 210 × 3.5 mm

Timer Options: 10 / 15 / 20 minutes

Material: ABS + Silicone



# innoGIO

Manufactured in China for **INNOGIO GLOBAL** sp. z o.o.  
Ostródzka street 74H, 03-289 Warsaw, Poland  
**[www.innogio.com](http://www.innogio.com)**



[@innoGIO.official](https://www.facebook.com/innoGIO.official)



[@innogio\\_com](https://www.instagram.com/innogio_com)

-  Scan the QR code to download the manual
-  Zeskanuj kod QR, aby pobrać instrukcję
-  Scannen Sie den QR-Code, um das Handbuch herunterzuladen
-  Scannez le code QR pour télécharger le manuel
-  Escanee el código QR para descargar el manual
-  Digitalize o código QR para descarregar o manual
-  Scansiona il codice QR per scaricare il manuale
-  Scan de QR-code om de handleiding te downloaden
-  Скануйте QR-код, щоб завантажити інструкцію
-  Naskenujte QR kód a stáhněte si návod k použití
-  Naskenujte QR kód a stiahnite si príručku
-  Scanează codul QR pentru a descărca manualul
-  Nuskaitykite QR kodą, kad atsisiųstumėte instrukciją
-  Nolasiet QR kodu, lai lejupielādētu rokasgrāmatu
-  Lae juhend alla, skannides QR-koodi
-  A kézikönyv letöltéséhez olvassa be a QR-kódot.
-  لي لڤل لى زنتل ءة ىرسلا ءباجتسالا زم ر حسم

